

## 01 Product Description

Product Name : Powder Free Sterile Latex Surgical Gloves for Gynaecological Use  
Material : Natural Rubber Latex  
Colour : Creamy White  
Shape : Anatomic  
Cuff : Beaded  
External Surface: Chlorinated  
Internal Surface : Chlorinated  
Size : 6.5(S), 7.5(M) and 8.5(L)  
Sterilization : Ethylene Oxide (EO) / GAMMA (R)  
Shelf Life : 3 years  
Basic UDI-DI : 806363LSGFM4

## 02 Intended Use

This disposal medical device is made up of natural rubber latex which is anatomically shaped with thumb position towards the Palm side of the index finger which reduces the fatigue on the hands, intended to be worn on the hands usually in surgical settings to provide barrier against potentially infectious fluids and other contaminants.

## 03 Product Classification & Standard compliance

Medical Device Classification : Class IIa, Rule#06  
Conformity assessment route : Annex II section 4 of council directive (EU) 2017/745  
Regulatory Authority : DNV Product Assurance As  
Notified Body Number : 2460  
Standard Compliance :  
EN 455-1:2020+A2:2024, EN 455-2:2024, EN 455-3:2023, EN 455-4:2009, EN ISO 15223-1:2021, EN 1041:2013, EN 566-1:2024, EN ISO 11135: 2014/A1:2019, EN ISO 11137-1:2015/A2:2019, EN ISO 11137-2:2015/A1:2023, EN ISO 11737-1:2018/A1:2021, EN ISO 11737-2:2020, ISO 11607-1:2020/A1:2023, ISO 11607-2:2020/A1:2023, EN ISO 10993-1:2020, EN ISO 10993-5:2009, EN ISO 10993-7:2008, EN ISO 10993-10:2021, EN ISO 10993-11:2017, EN ISO 10993-23:2021, ISO 10282:2023, ASTM D 3577-19(2023), ISO 13485:2016, ISO 14001:2015, ISO 9001 : 2015, IS 13422 : 2024.

## 04 Storage Instruction

Gloves must be stored in cool dry environment which is dust free. Cartons and Boxes must be stored unopened until required. Recommended storage temperature is 5°C-35°C. Avoid exposure to direct light, heat and excessive humidity. As ozone is deleterious, storage area should not contain any equipment capable of generating ozone such as ultraviolet light, fluorescent lights, mercury vapour lamps, photocopier, high voltage equipment, x-ray units, electric motors and electro surgical equipments.

## 05 Indication For Use

- ☞ Dry hands thoroughly before donning.
- ☞ Protective gloves should only be used for the intended application and in the correct size.
- ☞ These are sterile gloves for single use only.
- ☞ Users should take care when using the gloves. Using them solely according to their intended application.
- ☞ Before usage, inspect the gloves for any defect or imperfection.

## 06 Contraindication

- ☞ Latex gloves are made of Natural rubber latex, which may cause allergic reactions including anaphylaxis response if the user is allergic to latex.
- ☞ Gloves contain Natural Latex; persons who are sensitive to Latex should consult a physician before using.

## 07 Precautions


- ☞ Do not use if package is damaged or wet.
- ☞ Risk of reuse: May cause infection, allergic reaction and poor barrier protection.
- ☞ Gloves shall not be worn where there is a risk of entanglement by moving parts of machines is needed.
- ☞ The results do not reflect the actual duration of protection in the workplace due to other factors influencing the performance, such as temperature, abrasion, degradation etc.

## 08 Warnings

- ☞ Dispose off the devices and packaging after use as per Bio-Medical waste laws.
- ☞ Do not re-sterilize.
- ☞ The product contains Natural Rubber Latex which may cause allergic reactions including anaphylactic responses to some individuals.
- ☞ The gloves not intend to prevent Electrical shock care should be taken to have proper earthing in Medical Device Electrical appliances user.
- ☞ Necessary caution to be practiced against probable Electrical Hazards.


## 09 Directions for use

### Glove Opening and Donning Procedure :

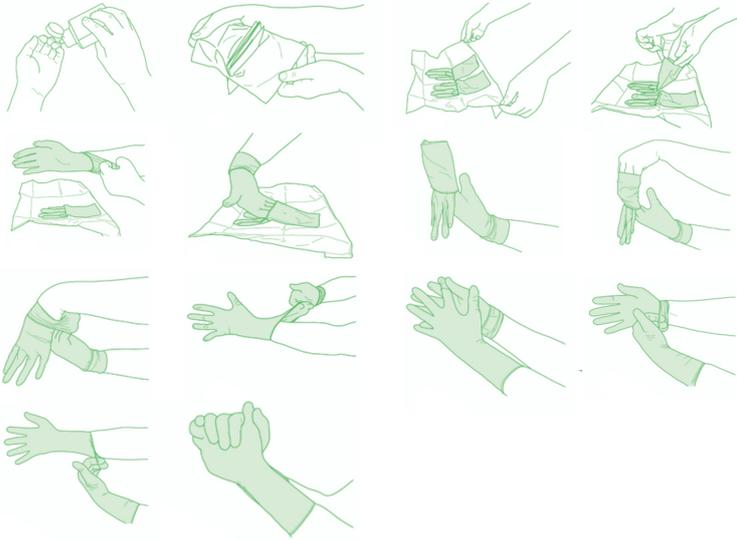

- ☞ Remove the walled gloves (inner wrapper) from the Pouch (outer wrapper) by Peel open from the corner for Paper Pouch (Peel down to open pouch).
- ☞ Open the walled glove to see “Left” and “Right” compartment.
- ☞ Pinch back upper and lower flaps of the inner wrapper.
- ☞ Using the middle flaps, open the wrapper touching only the 1 inch margin for safety.
- ☞ Be sure wrapper does not close over gloves after opening to avoid contamination.
- ☞ Using the thumb and the first two fingers of the non-dominant hand, pinch the cuff of the folded edge of the glove cuff for the dominant hand, touching only the inside surface of the glove.
- ☞ Slide dominant hand in to the gloves keeping hand point downwards and pull up to wrist.
- ☞ Using the glove hand insert the 4 fingers under the cuff of the other glove and pull the glove up to the arm.
- ☞ Adjust the gloves as necessary.

### Glove Removal Procedure :

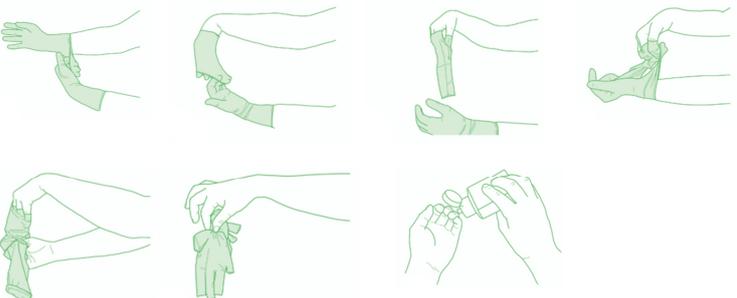

- ☞ Take hold of the first glove at the wrist.
- ☞ Fold it over and peel it back, turning it inside out as it goes. Once the glove is off, hold it with your gloved hand.
- ☞ To remove the other glove, place your bare fingers inside the cuff without touching the glove exterior. Peel the glove off from the inside, turning it inside out as it goes. Use it to envelope the other glove.

## 10 Explanation of Symbols

-  Manufacturer
-  Authorized representative in the European Community  
AMSTERMED B.V.,  
Saturnusstraat 46-62,  
Unit 032, 2132 HB Hoofddorp,  
The Netherlands
-  Date of Manufacture / Country of Manufacture
-  Use by date
-  Lot No
-  Reference Number / Catalogue Number
-  Serial Number
-  CE Logo
-  Sterilization using Ethylene oxide
-  Sterilization using Irradiation
-  Single Sterile Barrier System
-  Do not re-sterilize
-  Do not use package is damaged or wet
-  Keep away from sunlight
-  Keep dry
-  Temperature limit
-  This way up or end up
-  Keep away from Ozone
-  Single Use
-  Instruction for Use
-  Caution
-  Latex Caution
-  Medical Device
-  Unique Device Identifier

### 11 Personal Protective Equipment

EN ISO 374-1/Type A    EN ISO 374-1/Type B    EN ISO 374-1/Type C



UVWXYZ



XYZ



TYPE A Permeation min. 6 chemicals, level 2 (table 1)  
TYPE B Permeation min. 3 chemicals, level 2 (table 1)  
TYPE C Permeation min. 1 chemical, level 1 (table 1)

EN ISO 374-5:2016



Protective gloves protecting from viruses, bacteria and fungus

EN ISO 374-5:2016



Protective gloves protecting from bacteria and fungus

**Table 1** Performance Level

Chemicals	Time	Protection index
A Methanol	>10 min	1
B Acetone	>30 min	2
C Acetonitrile	>60 min	3
D Dichloromethane	>120 min	4
E Carbon disulphide	>240 min	5
F Toluene	>480 min	6
G Dimethylamine		
H Tetrahydrofuran		
I Ethyl acetate		
J n-Heptane		
K Sodium hydroxide 40%		
L Sulphuric acid 96%		
M Nitric acid 65%		
N Acetic acid 99%		
O Ammonium hydroxide 25%		
P Hydrogen peroxide 30%		
S Hydrofluoric acid 40%		
T Formaldehyde 37%		

EN ISO 374-1/Type B

EN ISO 374-5:2016



KLMP



VIRUS



**CE** In Compliance with the Harmonized Standards  
Personal Protective Equipment, risk category III under  
0598 Regulation (EU)2016/425

EN ISO 21420 : 2020  
EN ISO 374-1 : 2016  
EN ISO 374-2 : 2019  
EN 16523-1 : 2015+A1:2018  
EN ISO 374-4 : 2019  
EN ISO 374-5 : 2016

Regulatory Authority : SGS Fimko Oy,  
Takomotie 8, FI-00380, Helsinki, Finland.  
Notified Body Number : 0598

### 12 Warnings as per EN 374-1 & 5

- This information does not reflect the actual duration of protection in the work place and the differentiation between mixtures and pure chemicals.
- The Chemical resistance has been assessed under laboratory conditions from samples taken from the palm only (except in cases where the glove is equal to or over 400mm – where the cuff is tested also) and relates only to the chemical tested. It can be different if the chemical is used in a mixture.
- It is recommended to check that the gloves are suitable for the intended use because the conditions at the work place may differ from the type test depending on temperatures, abrasion and degradation.
- When used, protective gloves may provide less resistance to the dangerous chemical due to changes in physical properties.
- Movements, snagging, rubbing, degradation caused by the chemical contact etc., may reduce the actual use time significantly.
- For corrosive chemicals, degradation can be the most important factor to consider in selection of chemical resistant gloves.
- Before usage, inspect the gloves for any defect or imperfections.
- For Single use only.
- The Penetration resistance has been assessed under laboratory conditions and relates only to the tested specimen.

### 13 Test Results as per EN ISO 374-1 & 5

#### 1. EN ISO 21420 : 2020 Protective Gloves - General Requirement and test methods

400 mm				
Type	Declared Size	6.5	7.5	8.5
Glove	Length (mm)	406	396.3	401.7
	Circumference (mm)	169.7	189.3	208
	Dexterity	Performance Level 5		
	pH value	6.1		
480 mm				
Type	Declared Size	6.5(S)	7.5(M)	8.5(L)
Glove	Length (mm)	497	489	490
	Circumference (mm)	167	188.3	205
	Dexterity	Performance Level 5		
	pH value	6.0		

#### Polyaromatic Hydrocarbons content

	Results
Benz (a) anthracene	Not Detected
Chrysene	Not Detected
Benzo (b) Fluoranthene	Not Detected
Benzo (k) Fluoranthene	Not Detected
Benzo (a) pyrene	Not Detected
Dibenz (a,h) anthracene	Not Detected
Benzo (e) pyrene	Not Detected
Benzo (j) Fluoranthene	Not Detected

### 2. EN ISO 374-2 : 2019

Protective gloves against dangerous chemicals and micro-organisms - Part 2 : Determination of resistance to penetration

Resistance to Penetration	Status
Air Leakage	Pass
Water Leakage	Pass

Clause	Test Name	Test Results	Performance level
7.2	Air leak test (Air Pressure Used: 0.5 kPa)	Sample size	Pass
		Leakage	
		6.5 (S) 7.5 (M) 8.0 (L) 8.5 (L)	
		No Leakage No Leakage No Leakage No Leakage	
7.3	Water leak test	Sample size	Pass
		Leakage	
		6.5 (S) 7.5 (M) 8.0 (L) 8.5 (L)	
		No Leakage No Leakage No Leakage No Leakage	

### 3. EN ISO 374-4 : 2019

Protective gloves against dangerous chemicals and micro-organisms - Part 4 : Determination of resistance to degradation by chemicals

Chemical / CAS No	Exposure Duration	Test Results (Percentage change in puncture resistance)	Observation
Sodium hydroxide 40% 1310-73-2 PALM	60±5 minutes	Glove sample	No change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-4.5 -2.7 -6.3 -4.5 1.837	
Sodium hydroxide 40% 1310-73-2 CUFF	60±5 minutes	Glove sample	No change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-9.1 -8.9 -12.9 -10.3 1.1666	
Sulphuric acid 96% 7664-93-9 PALM	60±5 minutes	Glove sample	Severe Swelling & Colour change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		61.2 63.5 61.8 62.2 1.166	
Sulphuric acid 96% 7664-93-9 CUFF	60±5 minutes	Glove sample	Severe Swelling & Colour change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		58.7 60.7 57.0 58.8 1.828	
Nitric acid 65% 7697-37-2 PALM	60±5 minutes	Glove sample	Severe Swelling & Colour change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		54.0 53.0 52.8 53.3 0.630	
Nitric acid 65% 7697-37-2 CUFF	60±5 minutes	Glove sample	Severe Swelling & Colour change
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		59.5 57.6 57.3 58.1 1.187	
Acetic acid 99% 64-19-7 PALM	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		41.2 40.3 41.0 40.9 0.479	
Acetic acid 99% 64-19-7 CUFF	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		17.8 15.2 15.7 16.2 1.378	
Hydrogen peroxide 30% 7722-84-1 PALM	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-5.8 -4.2 -6.1 -5.3 0.989	
Hydrogen peroxide 30% 7722-84-1 CUFF	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-4.0 -3.6 -2.5 -3.4 0.753	
Formaldehyde 37% 50-00-0 PALM	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-5.4 -5.4 -2.7 -4.5 1.579	
Formaldehyde 37% 50-00-0 CUFF	60±5 minutes	Glove sample	Slight swelling
		Result (%)	
		1 2 3	
		Mean Standard Deviation	
		-10.7 -11.8 -8.7 -10.4 1.545	

### 4. EN 16523-1 : 2015+A1:2018

Determination of material resistance to permeation by chemicals - Part 1 : Permeation by Liquid chemical under conditions of continuous contact.

Chemical CAS No	Loop System /collection medium	Analytical Technique used	Mean thickness (mm)	NBT at NPR 1.0 µg cm <sup>-2</sup> min <sup>-1</sup> (minutes)	Performance level accordance to BS EN ISO 374-1:2016+A1:2018 Table 1	Observation
Sodium hydroxide 40% 1310-73-2 PALM	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.22	>480	6	No change
			0.21	>480		
			0.22	>480		
Sodium hydroxide 40% 1310-73-2 CUFF	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.16	>480	6	No change
			0.15	>480		
			0.16	>480		
Sulphuric acid 96% 7664-93-9 PALM	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.21	69	3	Severe swelling & Colour change
			0.22	67		
			0.21	65		
Sulphuric acid 96% 7664-93-9 CUFF	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.15	44	2	Severe swelling & Colour change
			0.16	40		
			0.15	42		
Nitric acid 65% 7697-37-2 PALM	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.22	136	4	Severe swelling & Colour change
			0.21	145		
			0.22	126		
Nitric acid 65% 7697-37-2 CUFF	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.15	62	3	Severe swelling & Colour change
			0.16	64		
			0.15	66		
Hydrogen peroxide 30% 7722-84-1 PALM	Closed loop /Grade 3 water	Periodic measurement with HPLC	0.22	340	5	Slight swelling
			0.21	346		
			0.22	348		
Hydrogen peroxide 30% 7722-84-1 CUFF	Closed loop /Grade 3 water	Periodic measurement with HPLC	0.15	300	5	Slight swelling
			0.16	305		
			0.15	310		
Formaldehyde 37% 50-00-0 PALM	Closed loop /Grade 3 water	Periodic measurement with HPLC	0.21	24	1	Slight swelling
			0.22	27		
			0.21	29		
Formaldehyde 37% 50-00-0 CUFF	Closed loop /Grade 3 water	Periodic measurement with HPLC	0.16	13	1	Slight swelling
			0.15	15		
			0.16	18		
Acetic acid 99% 64-19-7 PALM	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.21	6	0	Slight swelling
			0.22	7		
			0.21	6		
Acetic acid 99% 64-19-7 CUFF	Closed loop /Grade 3 water/350 rpm	Continuous measurement with conductivity electrode	0.15	3	0	Slight swelling
			0.16	4		
			0.15	3		

### 5. EN ISO 374-5:2016

Protective gloves against dangerous chemicals and microorganisms - Part 5 : Terminology and Performance Requirements for microorganisms risks

Protected against viruses according to ISO 16604 procedure B

Specimen	Thickness in mm	Observe the Presence of Plaques on Agar Plate (PFU)	Weight in gm	Results (Pass/Fail)	Assay Titer Ratio
1	0.21	<1	1.20	Pass	1.3
2	0.20	<1	1.22	Pass	
3	0.21	<1	0.19	Pass	
Negative control	No Plaque Formation observed				
Positive control	Plaque Formation observed				

Plate exposes during testing (Settle Plate): No plaque formation observed.

➤ Declaration of conformity is supplied with the product.